


MITRAL HICYCLE DURABILITY TESTER

Product Code: 27616

This enhanced version of the ViVidro HiCycle Durability Tester supports mitral valve durability testing for devices up to 55 mm in diameter. The Mitral HiCycle tests up to three devices simultaneously.

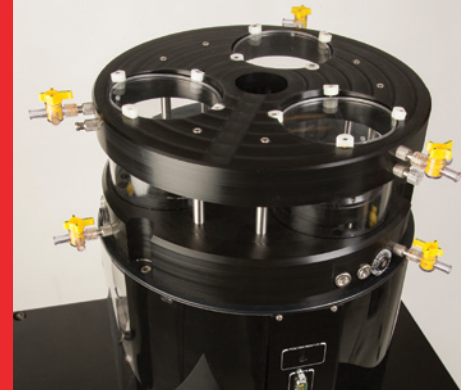

The Mitral HiCycle system uses the same electrodynamic motor as the original ViVidro HiCycle to drive three pistons within fluid-filled, temperature-controlled chambers moving the device rather than the fluid. The result is more stable fluid flow. Cycle rates are adjustable from 500 – 1,800 CPM (8.3 - 30 Hz).

The Mitral HiCycle is fully compatible with all the other ViVidro HiCycle components and accessories including:

- Data Acquisition Software
- Stroboscope
- HiCycle Cart
- Manometer
- Bellows

The Mitral HiCycle is available as a complete system or an upgrade kit can be purchased to convert existing HiCycles.



CONTACT US

HEAD OFFICE

ViVitro Labs Inc.
455 Boleskine Road
Victoria, BC, V8Z 1E7, Canada
1 250.388.3531
info@vivitrolabs.com
skype: Vivitro.labs.inc.
www.vivitrolabs.com

NORTH AMERICA

Douglas Bigrigg, B. Eng.
Major Accounts Manager
ViVitro Labs Inc.
455 Boleskine Road
Victoria, BC, V8Z 1E7, Canada
1 250.388.3531 x217
DBigrigg@vivitrolabs.com
www.vivitrolabs.com

LATIN AMERICA (EXCLUDING MEXICO)

Andreluiz de Sousa Cunha
Críons
R. Viamão 1267,
LJ8 Belo Horizonte, Brazil
30431-253 Brazil
+55(31) 2515-2505 / 2515-7017
andre@crions.com
www.crions.com

EUROPE, ISRAEL AND AFRICA

Karim Mouneimne
President & CEO
ProtomedLabs SASU
24-26 BD Gay Lussac
Immeuble Energie
F-13014 Marseille, France
+33 4 86 68 68 10
k.mouneimne@protomedlabs.com
www.protomedlabs.com

CHINA

Shining Y. W. Sun
沈阳福帝斯科技有限公司
Shenyang Fortis Tech. Ltd.
024-23286966
yuwei.sun@fortis-hitech.com
www.fortis-hitech.com

JAPAN

Toshiharu Miyake
CEO
Comfort Medical, Inc
Omagi 1751-1-101
Midori-ku, Saitama-shi
Saitama, Japan
+33 4 86 68 68 10
miyake@comfort-md.com
www.comfort-md.com

MIDDLE EAST

Mohan Das
Middle East Medical Supplies
PO Box 26143 Tahir Musbah Al Kendi Bldg.
Abu Dhabi, United Arab Emirates
+971 506221598
mohanda@memsuae.com
www.memsuae.com

INDIA

Harish Rao
Emergent Medi-Tech (I) Pvt. Ltd.
116-119 Damji Shamji Industrial Complex
9, L.B.S. Marg, Kurla (West)
Mumbai, 400070, India
91-98-20210150
91-22-67825555
info@emergent.in
www.emergent.in

No part of this publication may be reproduced, stored in a retrieval system or transmitted in any form or by any means, electronic, mechanical, photocopying, recording, scanning or otherwise without the prior written permission of the Publisher.

Limit of Liability/Disclaimer of Warranty: While the publisher and author have used their best efforts in preparing this document, they make no representation or warranties with respect to the accuracy or completeness of the contents of this document and specifically disclaim any implied warranties of merchantability or fitness for a particular purpose. No warranty may be created or extended by sales representatives or written sales material. The advice and strategies contained herein may not be suitable for your situation. Neither the Publisher nor author shall be liable for any loss of profit or any other commercial damages, including but not limited to special, incidental, consequential or other damages.

455 Boleskine Road, Victoria B.C. Canada, V8Z 1E7
info@ViVitroLabs.com | (250) 388-3531 | www.ViVitroLabs.com

Copyright © 2018 by ViVitro Labs Inc. All rights reserved.

 **ViVitro Labs Inc.**
First in Cardiovascular Device Testing