

# REAL-TIME WEAR TESTER (RWT)

Product Code: 20174

*The ViVitro Real-Time Wear Tester (RWT) is a quasi-real-time wear tester for heart valve replacements, venous Valve replacements, occluder devices, cardiac patches or strips, medical sutures, staples or clips, and other cardiovascular devices. It is customizable for devices and device mounting arrangements (including conduit mounted valves and large prosthetic heart valves).*

*The RWT gives users full control over the waveform driving the system while testing numerous devices simultaneously.*



With the Hitest Data Acquisition System, RWT operators are able to collect and analyze pressure data according to ISO5840 standards through an intuitive user interface. Quasi-real-time results may identify failure modes similar to animal tests for prosthetic heart valves.

## > Features



- Operator control over input waveform
- User interface designed for easy waveform download
- Controller stores up to 5 waveforms to drive the SuperPump
- Visible operation for easy viewing and image capture of testing in progress or real-time inspection
- Easy access cartridges allow quick valve change-out
- Tests up to six devices simultaneously
- Durability and fatigue testing up to 200 BPM

## > Testing Capabilities

- Heart Valve Replacements
- Venous Valve Replacements
- Occluder Devices
- Cardiac Patches or Strips
- Medical Sutures, Staples or Clips
- Other Cardiovascular Devices

## > Specifications

<b>Device Capacity</b>	1 - 6
<b>Dry Weight</b>	Controller: 8 kg Actuator + Tank: 46 kg
<b>Dimensions</b>	95 x 37 x 43 cm
<b>Tank Volume</b>	30 L
<b>Temperature Control</b>	up to 40°C ± 0.5°C
<b>Resettable Cycle Counter</b>	0 - 99,999 Million - resolution of 1,000 cycles
<b>Cycle Rate</b>	3 - 200 BPM

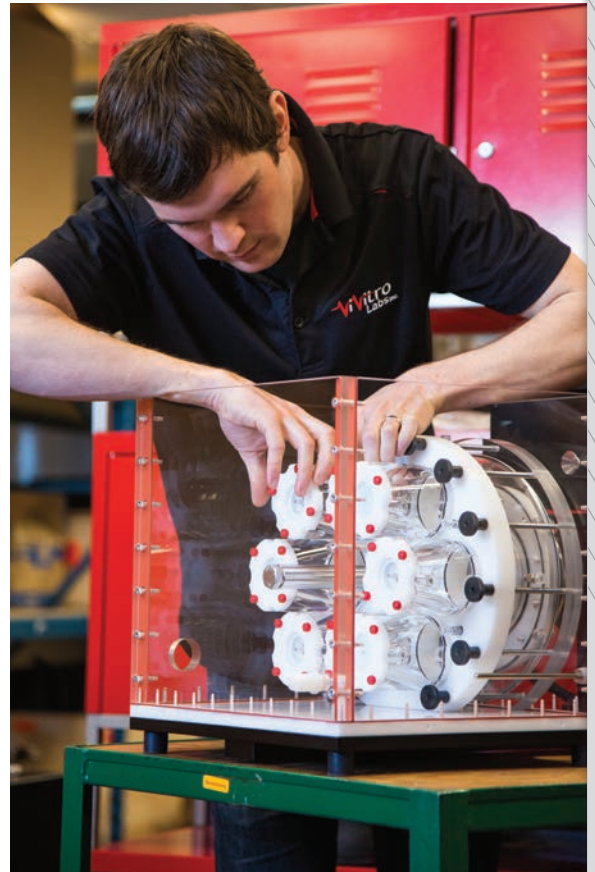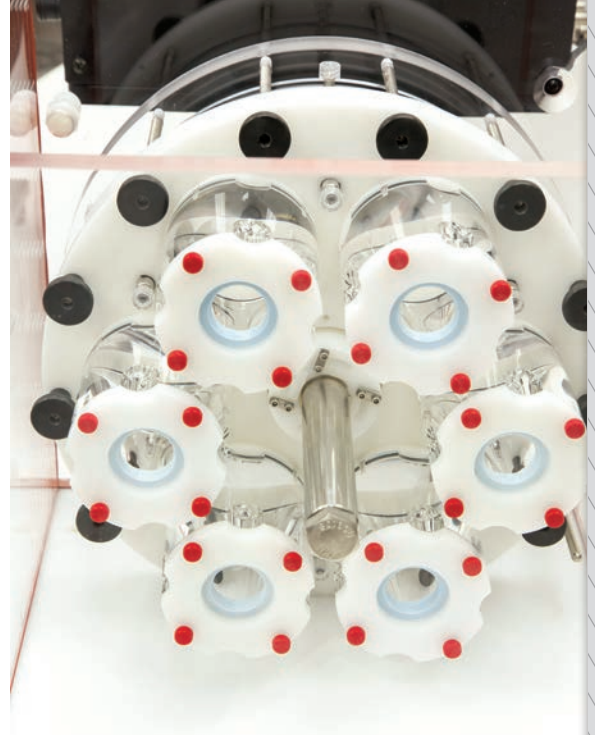

No part of this publication may be reproduced, stored in a retrieval system or transmitted in any form or by any means, electronic, mechanical, photocopying, recording, scanning or otherwise without the prior written permission of the Publisher.

**Limit of Liability/Disclaimer of Warranty:** While the publisher and author have used their best efforts in preparing this document, they make no representation or warranties with respect to the accuracy or completeness of the contents of this document and specifically disclaim any implied warranties of merchantability or fitness for a particular purpose. No warranty may be created or extended by sales representatives or written sales material. The advice and strategies contained herein may not be suitable for your situation. Neither the Publisher nor author shall be liable for any loss of profit or any other commercial damages, including but not limited to special, incidental, consequential or other damages.

455 Boleskine Road, Victoria B.C. Canada, V8Z 1E7 | (250) 388-3531  
24-26 Blvd. Gay Lussac, Immeuble Energie, 13014 Marseille, France | +33 4 86 68 68 10

info@ViVitroLabs.com | [www.ViVitroLabs.com](http://www.ViVitroLabs.com)

Copyright © 2023 by ViVitro Labs Inc. All rights reserved.

**Vivitro**  
Labs